

# AJEX SERIES

AJEX 160H

It's the multi-player

compact

powerful

versatile

cost effective



J-TYPE stand



M-TYPE stand




V-TYPE stand




AJEX Meditech, Ltd. since 1960

38-1, Mangpo-dong, Yeongtong-ku, Suwon-city, Kyunggido, Korea  
TEL : 82-31-204-4200~2 FAX : 82-31-204-4203  
email : ajextrad@kotis.net URL : www.ajex.co.kr

### What's difference

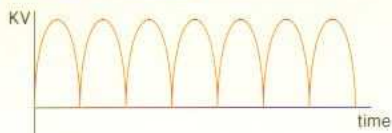

Conventional X-ray generator

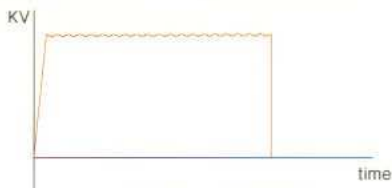

H F X-ray generator

### Features


- Applying high frequency inverter type generator (60kHz)
- Digital displays of kV, mAs and time
- Overload protection for kV, mA, Time and Temp
- Self-managible APR function with memory keys
- Soft touch user panel
- Wide range input power acceptable with electronic voltage compensator

### Advantages


- Powerful output covering most radiographic technique
- Superb quality images with fine focal spot X-ray tube
- Excellent reproducibility not affected by line voltage fluctuation
- Reduced the hazards of X-ray radiation
- Prolonged lifetime of X-ray tube
- Ordinary wall-outlet power use

### Specification


Items	Models	AJEX 160H
Max. output		3.3kW
Input Line Voltage		220V~, 1 Phase, 50/60Hz
Rating		40~100kV 0.6~135mAs (max. 60mA)
Exposure time		0.02~3.0sec
Focal spot		2.0mm
Memory keys		27(9x3)
B.L.D		100W(50W option) Collimators with rotatable adapter Full field light beam type collimators with 30sec. timer
Size(WXDXH)		253 x 341 x 248(mm) Exclude hand grip
Weight		17.5kgs
Bucky interface kit(option)		

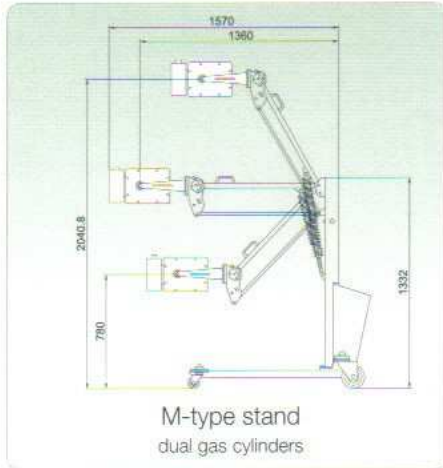

M-type stand  
dual gas cylinders

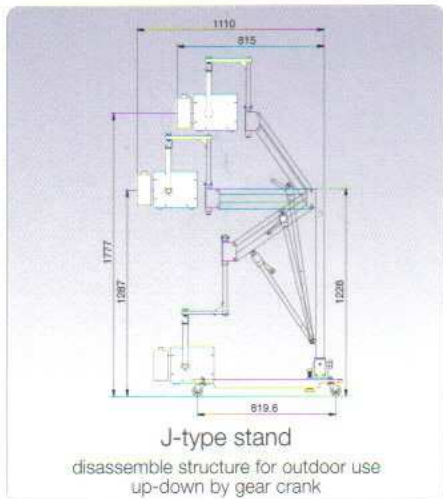

J-type stand  
disassemble structure for outdoor use  
up-down by gear crank

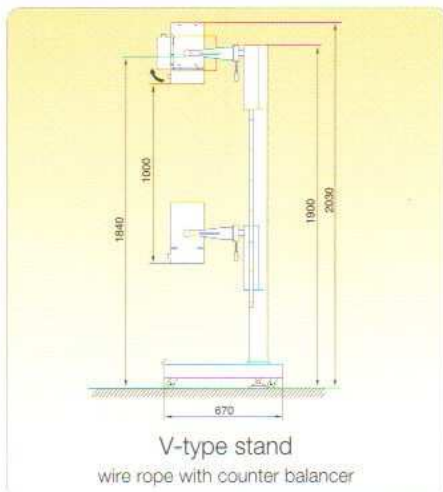

V-type stand  
wire rope with counter balancer

Distributed by

**Vega Group**  
Elektronik Tıbbi Cihazlar  
[www.vegagroup.com.tr](http://www.vegagroup.com.tr)

**AJEX** SERIES